

Quick Start Guide


Weight Scale

Your kit contains the following supplies:

- A.** Weight scale
- B.** 4 AAA batteries (preInstalled)

Getting started

- Make sure the Prolaio phone is turned on.
- The Home Screen will have an action card for a measurement task.
- You can also take measurement as needed.

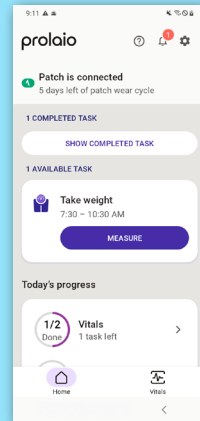


Weight measurement

- Place the scale on a flat, hard surface.
- Step onto the scale with your feet evenly planted on the surface.
- Wait a few moments for the scale to calculate your weight.




Weight measurement task

- Tap **Measure** on the measurement task action card.
- You will be presented with guided measurement instructions on your first measurement. Tap **Measure**.
- The next time, you can proceed directly to the measurement or use the guided instructions.
- Once your weight is calculated, you will be presented with a confirmation dialog. Tap **Confirm** to accept the measurement.

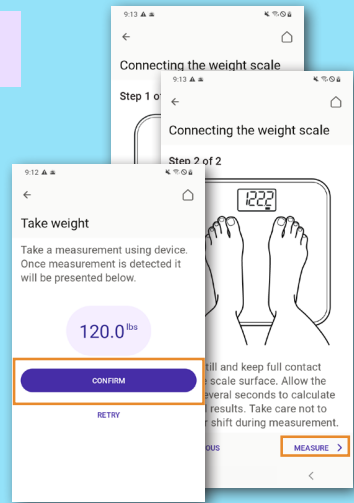

Weight measurement as needed

- Tap the **Weight Scale** action card from the Home Screen.
- Tap **Take Measurement**.
- Follow the same steps as described above.

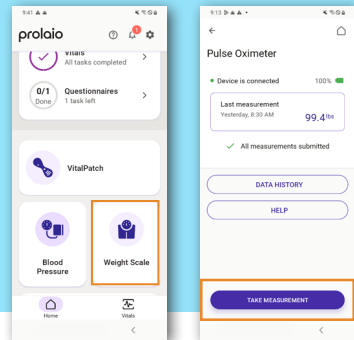

Common Error Codes

Error	Description	Solution
	Overload / too much weight is on the scale. The unit will shut off.	This scale is not suitable for weight in excess of 203kg/447.6lbs; cease using this scale.
	Low battery	Replace all four (4) AAA batteries with new batteries.
	Data upload has failed	Ensure your Prolaio phone is on and near the weight scale. Attempt to sync your results again.

